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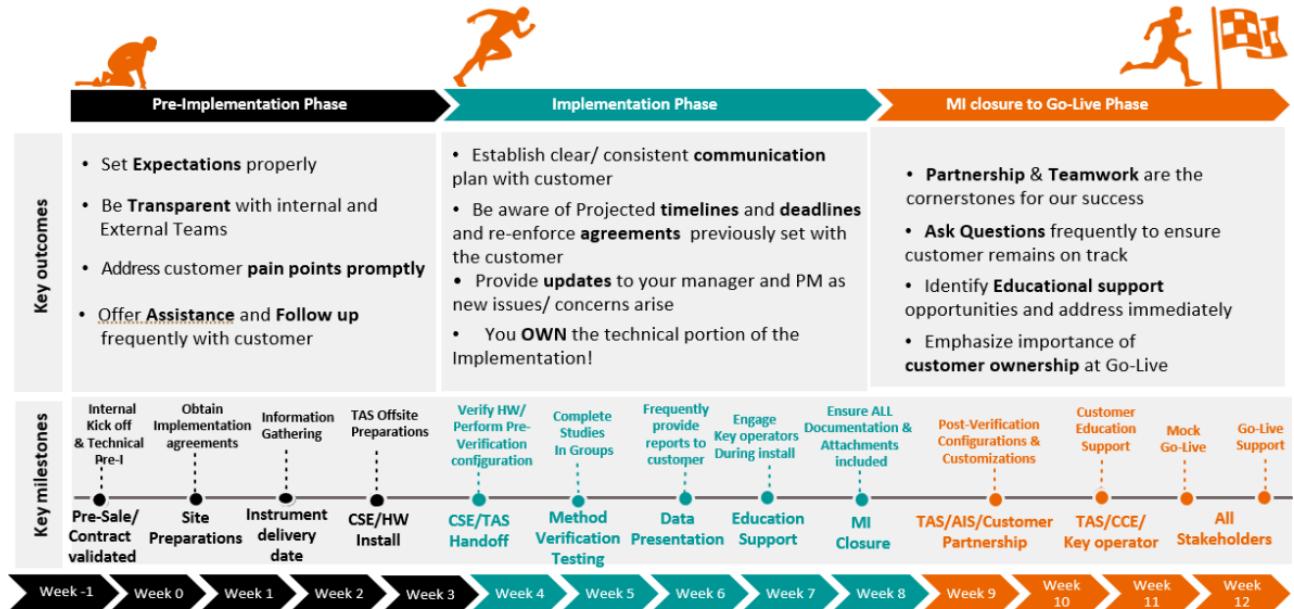
 SH1-SCI1

 LabManager

Select Tests
Assay Identification
Advanced

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Setup > Test Definition > IM Test Definition

Select Tests	Reference	Delta Check	Check	Concentration	Dilution Point	Linearity	Interpretation																																																																																																						
AFP (ng/ml)																																																																																																													
<input type="radio"/> Definition																																																																																																													
<input type="radio"/> Calculation																																																																																																													
<input type="radio"/> Calibration																																																																																																													
<input type="radio"/> Reagents																																																																																																													
<input type="radio"/> Lots																																																																																																													
<input checked="" type="radio"/> Ranges*																																																																																																													
<input type="radio"/> Repeat																																																																																																													
<table border="1"> <tr> <td>Range Name RR</td> <td colspan="7"> <input type="button" value="Add Range"/> <input type="button" value="Edit Range"/> <input type="button" value="Copy Range"/> <input type="button" value="Delete Range"/> </td> </tr> <tr> <td colspan="2">Defined Ranges</td> <td colspan="6"> Details <table border="1"> <tr> <td>Limits</td> <td>Critical Low 1.0000</td> <td>Low < 4.0000</td> <td>Normal 4.0000 - 5.0000</td> <td>High > 5.0000</td> <td>Critical High 6.0000</td> </tr> <tr> <td>Flag result</td> <td><input type="checkbox"/></td> <td><input checked="" type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input checked="" type="checkbox"/></td> </tr> <tr> <td>Repeat</td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> </tr> <tr> <td colspan="6">Number of Repeats</td> </tr> <tr> <td>Dilution</td> <td>Original</td> <td colspan="5"></td> </tr> <tr> <td>Replicates</td> <td>6</td> <td colspan="5"></td> </tr> <tr> <td>Reflex</td> <td>ALT</td> <td colspan="5"></td> </tr> <tr> <td>Dilution</td> <td>x1(Undiluted)</td> <td colspan="5"></td> </tr> <tr> <td>Replicates</td> <td></td> <td colspan="5"></td> </tr> <tr> <td>Comment</td> <td>Critical High</td> <td colspan="5"></td> </tr> </table> </td> </tr> <tr> <td colspan="8"> Demographics Information <table border="1"> <tr> <td>Specimen All</td> <td>Location All</td> <td>Age (Days)</td> <td>Sex</td> </tr> <tr> <td>Species All</td> <td>Status All</td> <td><input checked="" type="checkbox"/> All <input type="checkbox"/> From <input type="checkbox"/> To</td> <td><input type="checkbox"/> All <input type="checkbox"/> Not Specified <input checked="" type="checkbox"/> Male <input type="checkbox"/> Female</td> </tr> <tr> <td></td> <td></td> <td>Years</td> <td></td> </tr> </table> </td> </tr> </table>								Range Name RR	<input type="button" value="Add Range"/> <input type="button" value="Edit Range"/> <input type="button" value="Copy Range"/> <input type="button" value="Delete Range"/>							Defined Ranges		Details <table border="1"> <tr> <td>Limits</td> <td>Critical Low 1.0000</td> <td>Low < 4.0000</td> <td>Normal 4.0000 - 5.0000</td> <td>High > 5.0000</td> <td>Critical High 6.0000</td> </tr> <tr> <td>Flag result</td> <td><input type="checkbox"/></td> <td><input checked="" type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input checked="" type="checkbox"/></td> </tr> <tr> <td>Repeat</td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> </tr> <tr> <td colspan="6">Number of Repeats</td> </tr> <tr> <td>Dilution</td> <td>Original</td> <td colspan="5"></td> </tr> <tr> <td>Replicates</td> <td>6</td> <td colspan="5"></td> </tr> <tr> <td>Reflex</td> <td>ALT</td> <td colspan="5"></td> </tr> <tr> <td>Dilution</td> <td>x1(Undiluted)</td> <td colspan="5"></td> </tr> <tr> <td>Replicates</td> <td></td> <td colspan="5"></td> </tr> <tr> <td>Comment</td> <td>Critical High</td> <td colspan="5"></td> </tr> </table>						Limits	Critical Low 1.0000	Low < 4.0000	Normal 4.0000 - 5.0000	High > 5.0000	Critical High 6.0000	Flag result	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	Repeat	<input type="checkbox"/>	Number of Repeats						Dilution	Original						Replicates	6						Reflex	ALT						Dilution	x1(Undiluted)						Replicates							Comment	Critical High						Demographics Information <table border="1"> <tr> <td>Specimen All</td> <td>Location All</td> <td>Age (Days)</td> <td>Sex</td> </tr> <tr> <td>Species All</td> <td>Status All</td> <td><input checked="" type="checkbox"/> All <input type="checkbox"/> From <input type="checkbox"/> To</td> <td><input type="checkbox"/> All <input type="checkbox"/> Not Specified <input checked="" type="checkbox"/> Male <input type="checkbox"/> Female</td> </tr> <tr> <td></td> <td></td> <td>Years</td> <td></td> </tr> </table>								Specimen All	Location All	Age (Days)	Sex	Species All	Status All	<input checked="" type="checkbox"/> All <input type="checkbox"/> From <input type="checkbox"/> To	<input type="checkbox"/> All <input type="checkbox"/> Not Specified <input checked="" type="checkbox"/> Male <input type="checkbox"/> Female			Years					
Range Name RR	<input type="button" value="Add Range"/> <input type="button" value="Edit Range"/> <input type="button" value="Copy Range"/> <input type="button" value="Delete Range"/>																																																																																																												
Defined Ranges		Details <table border="1"> <tr> <td>Limits</td> <td>Critical Low 1.0000</td> <td>Low < 4.0000</td> <td>Normal 4.0000 - 5.0000</td> <td>High > 5.0000</td> <td>Critical High 6.0000</td> </tr> <tr> <td>Flag result</td> <td><input type="checkbox"/></td> <td><input checked="" type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input checked="" type="checkbox"/></td> </tr> <tr> <td>Repeat</td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> </tr> <tr> <td colspan="6">Number of Repeats</td> </tr> <tr> <td>Dilution</td> <td>Original</td> <td colspan="5"></td> </tr> <tr> <td>Replicates</td> <td>6</td> <td colspan="5"></td> </tr> <tr> <td>Reflex</td> <td>ALT</td> <td colspan="5"></td> </tr> <tr> <td>Dilution</td> <td>x1(Undiluted)</td> <td colspan="5"></td> </tr> <tr> <td>Replicates</td> <td></td> <td colspan="5"></td> </tr> <tr> <td>Comment</td> <td>Critical High</td> <td colspan="5"></td> </tr> </table>						Limits	Critical Low 1.0000	Low < 4.0000	Normal 4.0000 - 5.0000	High > 5.0000	Critical High 6.0000	Flag result	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	Repeat	<input type="checkbox"/>	Number of Repeats						Dilution	Original						Replicates	6						Reflex	ALT						Dilution	x1(Undiluted)						Replicates							Comment	Critical High																																													
Limits	Critical Low 1.0000	Low < 4.0000	Normal 4.0000 - 5.0000	High > 5.0000	Critical High 6.0000																																																																																																								
Flag result	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>																																																																																																								
Repeat	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>																																																																																																								
Number of Repeats																																																																																																													
Dilution	Original																																																																																																												
Replicates	6																																																																																																												
Reflex	ALT																																																																																																												
Dilution	x1(Undiluted)																																																																																																												
Replicates																																																																																																													
Comment	Critical High																																																																																																												
Demographics Information <table border="1"> <tr> <td>Specimen All</td> <td>Location All</td> <td>Age (Days)</td> <td>Sex</td> </tr> <tr> <td>Species All</td> <td>Status All</td> <td><input checked="" type="checkbox"/> All <input type="checkbox"/> From <input type="checkbox"/> To</td> <td><input type="checkbox"/> All <input type="checkbox"/> Not Specified <input checked="" type="checkbox"/> Male <input type="checkbox"/> Female</td> </tr> <tr> <td></td> <td></td> <td>Years</td> <td></td> </tr> </table>								Specimen All	Location All	Age (Days)	Sex	Species All	Status All	<input checked="" type="checkbox"/> All <input type="checkbox"/> From <input type="checkbox"/> To	<input type="checkbox"/> All <input type="checkbox"/> Not Specified <input checked="" type="checkbox"/> Male <input type="checkbox"/> Female			Years																																																																																											
Specimen All	Location All	Age (Days)	Sex																																																																																																										
Species All	Status All	<input checked="" type="checkbox"/> All <input type="checkbox"/> From <input type="checkbox"/> To	<input type="checkbox"/> All <input type="checkbox"/> Not Specified <input checked="" type="checkbox"/> Male <input type="checkbox"/> Female																																																																																																										
		Years																																																																																																											

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Definition / General / Advanced

The screenshot shows the 'CH Test Definition' configuration window divided into three main sections: 'General', 'Advanced', and 'Below Measuring Interval Configuration'.

General Tab:

- Select Test:** GluH_3 (mp/dL)
- Definition:** Definition*
- Calculation:**
- Calibration:**
- Mitigation:**
- Ranges:**
- Repeat:**
- Test Name / ID:** GluH_3
- Display Name:** GluH_3
- Print Name:** GluH_3
- Test Type:** IMT
- LIS Code:** GluH_3
- Test Version:** CH CI CH-HV
- LOINC:** 1.0 1.1 1.0
- Result Review Mode:** None
- Reuse Result for Ratio:**
- Result Time Limit (hours):** 99
- Analyte Stability (minutes):** 0

Advanced Tab:

- Troubleshooting:** Kinetic Mode
- Measuring Intervals:**
 - Serum:** Low: 50.00, High: 200.00
 - Plasma:** Low: 50.00, High: 200.00
 - Urine:** Low: 50.00, High: 200.00
- Repeat when Outside Measuring Interval:**
- Above Measuring Interval Configuration:**

Sample Volume: 50.00 uL	Dilution Factor: 4.00	Total Volume: 200.00 uL
-------------------------	-----------------------	-------------------------
- Below Measuring Interval Configuration:**

Sample Volume: 0 uL	Dilution Factor: 0	Total Volume: 0 uL
---------------------	--------------------	--------------------

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Select Tests GluH_3 (mg/dL)

- Definition
- Calculation
- Calibration
- Mitigation
- Ranges***
- Repeat

Reference	Delta Check	Check	Concentration	Dilution Point	Linearity	Interpretation																																								
Range Name																																														
Defined Ranges																																														
Details <table border="1"> <tr> <td>Limits</td> <td>Below Range</td> <td>Within Range</td> <td>Above Range</td> </tr> <tr> <td>Flag result</td> <td><input checked="" type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input checked="" type="checkbox"/></td> </tr> <tr> <td>Repeat</td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> </tr> <tr> <td>Number of Repeats</td> <td colspan="3"></td> </tr> <tr> <td>Dilution</td> <td>AUTO ...</td> <td colspan="2">AUTO +</td> </tr> <tr> <td>Replicates</td> <td colspan="3"></td> </tr> <tr> <td>Reflex</td> <td colspan="3"></td> </tr> <tr> <td>Dilution</td> <td colspan="3"></td> </tr> <tr> <td>Replicates</td> <td colspan="3"></td> </tr> <tr> <td>Comment</td> <td colspan="3"></td> </tr> </table>							Limits	Below Range	Within Range	Above Range	Flag result	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	Repeat	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Number of Repeats				Dilution	AUTO ...	AUTO +		Replicates				Reflex				Dilution				Replicates				Comment			
Limits	Below Range	Within Range	Above Range																																											
Flag result	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>																																											
Repeat	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>																																											
Number of Repeats																																														
Dilution	AUTO ...	AUTO +																																												
Replicates																																														
Reflex																																														
Dilution																																														
Replicates																																														
Comment																																														
<input type="button" value="Add Range"/> <input type="button" value="Edit Range"/> <input type="button" value="Copy Range"/> <input type="button" value="Delete Range"/>																																														
<input type="button" value="Reset Ranges"/> <input type="button" value="Save"/> <input type="button" value="Cancel"/>																																														

Select Tests GluH_3 (mg/dL)

- Definition
- Calculation
- Calibration
- Mitigation
- Ranges***
- Repeat

Reference	Delta Check	Check	Concentration	Dilution Point	Linearity	Interpretation																																								
Range Name																																														
Defined Ranges																																														
Details <table border="1"> <tr> <td>Limits</td> <td>Overdilution Point</td> <td>Within Range</td> <td>Dilution Point</td> </tr> <tr> <td>Flag result</td> <td><input checked="" type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> </tr> <tr> <td>Repeat</td> <td><input checked="" type="checkbox"/></td> <td><input type="checkbox"/></td> <td><input type="checkbox"/></td> </tr> <tr> <td>Number of Repeats</td> <td colspan="3"></td> </tr> <tr> <td>Dilution</td> <td>AUTO ...</td> <td colspan="2">AUTO +</td> </tr> <tr> <td>Replicates</td> <td>4</td> <td colspan="2"></td> </tr> <tr> <td>Reflex</td> <td colspan="3"></td> </tr> <tr> <td>Dilution</td> <td colspan="3"></td> </tr> <tr> <td>Replicates</td> <td colspan="3"></td> </tr> <tr> <td>Comment</td> <td colspan="3"></td> </tr> </table>							Limits	Overdilution Point	Within Range	Dilution Point	Flag result	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Repeat	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	Number of Repeats				Dilution	AUTO ...	AUTO +		Replicates	4			Reflex				Dilution				Replicates				Comment			
Limits	Overdilution Point	Within Range	Dilution Point																																											
Flag result	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>																																											
Repeat	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>																																											
Number of Repeats																																														
Dilution	AUTO ...	AUTO +																																												
Replicates	4																																													
Reflex																																														
Dilution																																														
Replicates																																														
Comment																																														
<input type="button" value="Add Range"/> <input type="button" value="Edit Range"/> <input type="button" value="Copy Range"/> <input type="button" value="Delete Range"/>																																														
<input type="button" value="Reset Ranges"/> <input type="button" value="Save"/> <input type="button" value="Cancel"/>																																														

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Select Tests
GluH_3 (mg/dL)

- Definition
- Calculation
- Calibration
- Mitigation
- Ranges
- Repeat*

Repeat Conditions for System Flags				
Repeat	Flags	Dilution	Replicates	
<input type="checkbox"/>	Clog	Original	1	
<input type="checkbox"/>	Sample Integrity	Original	1	
<input type="checkbox"/>	Insufficient Reagent	Original	1	
<input type="checkbox"/>	Abnormal Reaction	Original	1	
<input type="checkbox"/>	Measurement Error	Original	1	
<input type="checkbox"/>	Insufficient Sample	Original	1	

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Setup > SH & Sorting Configuration

To make configuration changes, stop all Sample Handlers

Error Sort Groups Selected Sort Group : Lab Default Group

Lab Default Group Add New Sort Group

HIL Error Group

Add New Sort Group +

Errors detected by Tube Vision System

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Install Defaults (кратко):

- USB Options – Enabled
- LIS Connection – Disabled
- LAS Connection – Disabled
- Patient Results Review Mode – Review None
- Send Results to QC Statistics – Enabled
- Enable Autorepeat – Disabled
- Automatic Results Reports – Enabled (Send to File → F:\UserRoot\Print To File Reports)
- Automatic Calibration Results Reports – Enabled (Send to File → тот же путь)

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- **HIL** – Manual HIL → Save

Calibration Setup:


- Mode: By System
- Automatically Order Calibrations: Disabled
- Perform QC with Calibration: Enabled
- Calibration Results: Automatic Acceptance Enabled
- Calibration Due: 1 hour → Save

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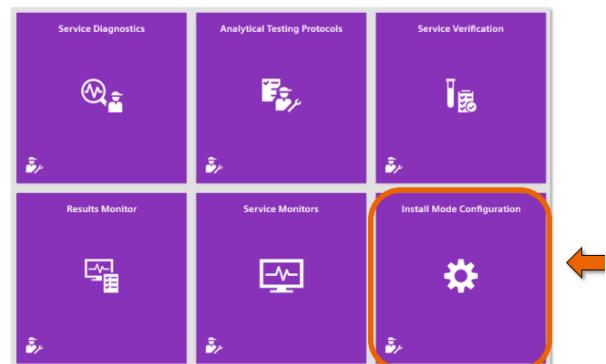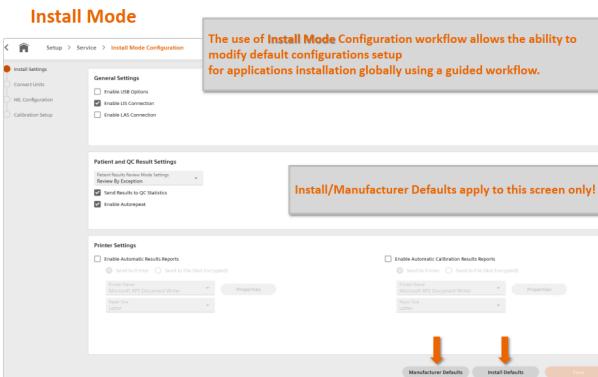
- Enable USB Options
- Enable LIS Connection
- Enable LAS Connection

Patient and QC Result Settings

Patient Results Review Mode Settings: Review None

- Send Results to QC Statistics
- Enable Autorepeat

Manufacturer Defaults

General Settings

- Enable USB Options
- Enable LIS Connection
- Enable LAS Connection

Patient and QC Result Settings

Patient Results Review Mode Settings: Review By Exception

- Send Results to QC Statistics
- Enable Autorepeat

Install Settings

Install Defaults

Printer Settings

- Enable Automatic Results Reports
 - Send to Printer
 - Send to File (Not Encrypted)

Manufacturer Defaults

Printer Settings

- Enable Automatic Results Reports
 - Send to Printer
 - Send to File (Not Encrypted)

Install Mode: Convert Units – Confirmation

Install Mode: Convert Units – Log

SIEMENS
Healthineers

HIL Configuration

Deactivate Install Mode

Calibration Setup

Calibration Order Settings

- Mode: By System
- Calibration Order Calibration (Enabled)
- Order QC with Calibration (Enabled)
- Auto Calibration Time: 100 ms (200-400)

Calibration Results

- Mode: By Assay Type
- QC Results Must Be Within Range (Enabled)
- QC Results Accuracy (Enabled)

Calibration Order Settings

- Calibration Order: 1

Calibration

- Cal-QC Storage Inventory
- Calibration Orders
- Calibrations Due
- Calibration Results
- Calibrator Alarms
- Calibrator Needs
- HIL Calibration

RECOMMENDATIONS:
Select Disable Perform QC with Calibration and then Save
Upon completion of Application Installation Configuration , settings should be changed according to the laboratory needs.
Changes can be done within this workflow or as part of the Calibration Setup

As per Atellica Solution | Assay Verification and Optional Configuration
Print No. LDAT-060.815.11.01.02, Revision date: 12.22

Deactivate Install Mode

- Must review Installation data first.
 - Then change setting from **Install Defaults** to **Manufacturer Default**.
 - You must be logged in as CSE to change the settings.
- Ensure all modules are in Stop mode to activate settings.
 - Navigate to Service > Install Mode Configuration.
 - At the bottom of the screen, click **Manufacturer Defaults**.
 - At the Confirmation window, click **Yes**.
 - Confirm all settings with the customer and customize as needed.
 - Click **Save**.

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Ручной ввод (Alliance Applications, Drugs of Abuse, User-Defined Assays):

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Service > Service Access > Service Verification > Assay Verification

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
Service Access > System Verification > Assay Verification Studies

SIEMENS
Healthineers

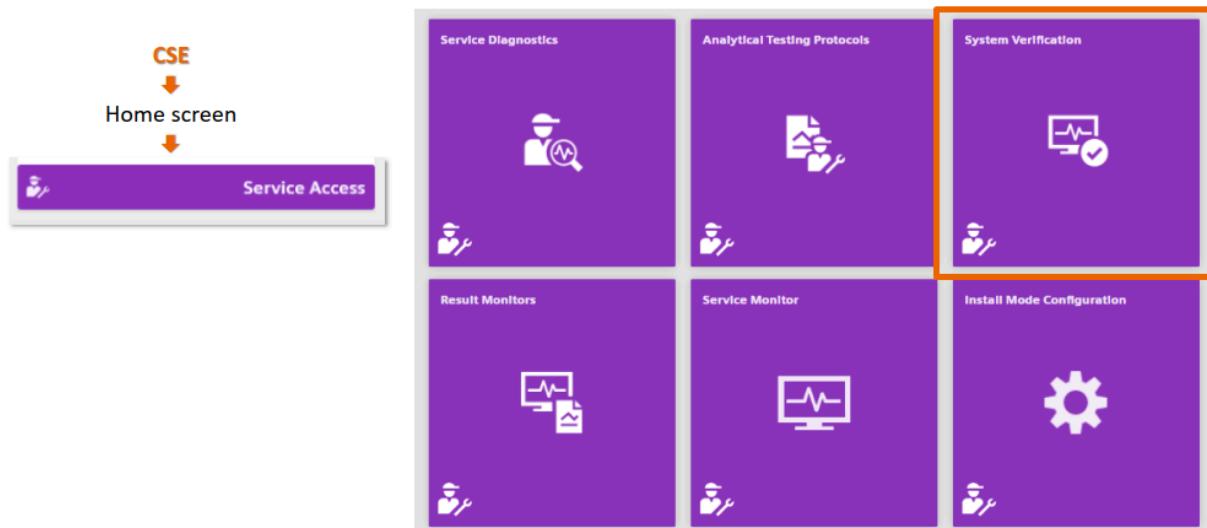Service Verification > Assay Verification

SIEMENS
Healthineers

The screenshot shows the 'System Verification Evaluations' section. It includes two radio button options: 'SVK Run A' and 'SVK Run B'. Below these are sections for 'Assay Verification Studies' and 'Assay Verification'. A dashed orange box highlights the 'Assay Verification Studies' section, which contains a radio button for 'Assay Verification' and a detailed description of its purpose. At the bottom right, there are 'Next' and 'Cancel' buttons.

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Test Definition > Operator-defined.

The screenshot shows the 'Operator-defined' test definition setup in the Atellica Solution Online Help. The interface is divided into several sections:

- Left Sidebar:** Shows a tree view of test definitions, with 'Definition' selected. Other options like 'Sample/Reagent', 'Calculation', 'High-dose Hook Effect', 'Calibration', 'Mitigation', 'Ranges', and 'Repeat' are also listed.
- Top Bar:** Includes icons for home, setup, test definition, and operator-defined, along with a search bar and user information.
- Central Area:**
 - Assay Identification:** Fields for Assay Name*, Units*, ID, Test version*, Display Name*, Print Name*, LIS Code*, LOINC, Investigation Only, and Kinetic Mode.
 - Measuring Intervals:** A table for Serum, Plasma, Urine, CSF, and RBC Hemolysate. For Serum, 'Low' is set to 4 and 'High' to 700. There is also a checkbox for 'Repeat when Outside Measuring Interval'.
 - Sample:** Predilution section with Sample Volume (µL) set to 50, Dilution Factor to 5, and Total Volume to 250µL.
 - Below Measuring Interval Configuration:** Fields for Sample Volume (µL), Dilution Factor, and Total Volume (µL).
 - Above Measuring Interval Configuration:** Fields for Sample Volume (µL), Dilution Factor, and Total Volume (µL).
 - General:** Result Review Mode (set to Now), Reuse Result for Ratio, Result Time Limit (Hours) set to 99, and Analyte Stability (minutes) set to 0.
- Bottom Buttons:** Create New, Create from Existing, Results Settings, Delete, Save (highlighted in orange), and Cancel.

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Test Definition > Off-System / Panels / Ratios.

The screenshot shows a list of existing tests on the left and a configuration dialog for a new test on the right. The new test is named 'GluH42' with LIS code 'GluH42'. It has a calculation rule where the low limit is marked as required (*). The enable/disable status is set to 'Enable'.

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The screenshot shows a list of tests with their current positions. The test 'aHAVM' is currently selected (indicated by a checked checkbox). The 'Move Selection' button is active, showing a range from -2 to 2. The 'Save' and 'Cancel' buttons are at the bottom.

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При необходимости заполните: - Physicians, - Locations, - Patient Status, - Species, - Sample Comment, - Patient Comment, - Result Comment.

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Перейдите: **Navigator > General Setup > Calibration Setup**.

Calibration Ordering

- Ordering Mode: By System
- Automatic Order Calibrations: Enabled

Order QC with Calibration

- Perform QC with Calibration: Disabled

Calibration Time

- Time: 01:00 AM

Calibration Results

- Results Mode: By System
- Automatic Acceptance: Enabled

QC Results Must be Within Ranges

Calibration Due Configuration

- Calibration Due Warning (hours): 1

Buttons: Restore Default, Save, Cancel

Calibration Ordering

- Ordering Mode: By Calibration Definition
- Calibrators Configured to Store
- Automatic Order Calibrations: Enabled
- Calibrators Configured Not to Store
- Automatic Order Calibrations: Disabled

Order QC with Calibration

- Perform QC with Calibration: Disabled

Calibration Time

- Time: 01:00 AM

Calibration Results

- Results Mode: By Assay Type
- CH
- Automatic Acceptance: Enabled
- QC Results Must be Within Ranges
- IM
- Automatic Acceptance: Enabled
- QC Results Must be Within Ranges

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The screenshot shows the QC Settings page in the LabManager software. It includes the following sections:

- Send Results to QC Statistics:** Options for sending results to QC or Patient.
- QC Mode Settings:** Selection between By System or By Assay Type. Sub-options include Never Disable, Enable, and Disable.
- Patient QC Flagging:** Option to Suspend Patient QC Flagging.
- RIIIBÄK Feature:** Selection between Off and On.
- QC Results Review Mode Settings:** Mode selection: Review By Exception.
- Hold Results by Analyzer:** List of analyzers (CH-01, CH-02, CH-03, CI-01, IM-01) with checkboxes for holding results.
- QC Schedule:** Option to Warn prior to QC Schedule if onboard control is unavailable.
- When to Warn:** Scheduling options for QC by Day and Time (Hours: 2, Minutes: 30) and QC Scheduled by Test Count (% Remaining: 25).

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The screenshot shows the 'QC Rules' tab selected under 'Setup'. It includes sections for 'Single Point Rules' and 'Multipoint Rules' with various configuration options like 'Error', 'Warning', 'SDa', and 'Off'. There are also sections for 'Custom single-point rule' and 'Custom multipoint rule'. At the bottom, there are two radio buttons: 'Evaluate regardless of a 1-2s rule error' (selected) and 'Evaluate only after a 1-2s rule error'. Below these are 'Restore Defaults', 'Refresh', 'Save', 'Reset', and 'Close' buttons.

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The screenshot shows the QC Scheduler interface. At the top, there are tabs for 'By Day and Time', 'By Test Count', and 'By Control'. Below the tabs, there's a search bar and a date selector showing 'Next Scheduled QC: 8:00 AM'. On the right side, there are checkboxes for 'Show Maintenance Schedule' and 'Show Calibration Schedule'. The main area is a grid where each cell represents a specific time slot on a specific day. The grid contains various tasks such as QC, Maintenance, and Calibration. A modal window is open for Sunday, titled 'View / Edit Schedule'. It shows a dropdown for 'Day' set to 'Monday' and 'Time' set to '08 : 00 AM'. There are two radio buttons: 'By Control Material' (selected) and 'By Panel'. A dropdown for 'Select Analyzer' is set to 'All'. Below these are tables for 'Material Name', 'Levels', and 'Tests' for various analyzers like ALPQUAL, CO2QUAL, CHEMQUAL, etc. At the bottom of the modal are 'Save', 'Cancel', and 'Copy Day' buttons.

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The screenshot shows the QC Scheduler interface with the 'By Test Count' tab selected. The main table lists tests (Test, Module, Test Count, Definitions, Operator Reset, Last Order, State) and actions (View/Edit, Restore Counter, Disable, Delete). A modal window titled 'View / Edit Schedule - By Test Count' is open for the 'Test Settings' of the 'AFP' test under 'Module CH-01'. The modal contains fields for 'Material Name' (DRUGQUAL), 'Levels' (1,2,3,4,5), and 'Assay' (AFP, CEA, Dig, eE2, hCG, TSH3-UL). Buttons for 'Save' and 'Cancel' are at the bottom right of the modal.

Test	Module	Test Count	Definitions	Operator Reset	Last Order	State
AFP	CH-01	0 <div style="width: 150px;"></div> 150	DRUGQUAL			Enabled
ALB	CH-02	0 <div style="width: 300px;"></div> 300	MULTIQUAL			Enabled
AFP	IM-01	0 <div style="width: 150px;"></div> 150	DRUGQUAL			Enabled
CAL	IM-02	0 <div style="width: 300px;"></div> 300	CHEMQUAL			Disabled
Carb	System	0 <div style="width: 300px;"></div> 300	ENZQUAL			Enabled

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The screenshot shows the QC Scheduler interface with the 'By Control' tab selected. The main table lists control materials (Control Name, Control Lot, Control ID, Control Lot Expiration, Assay, Usage, Definitions, Status, Last Order Submitted) and actions (Disable, Delete). A modal window titled 'Edit Control Lot' is open for the 'CHEMQUAL' control material. The modal contains fields for 'Control Name' (CHEMQUAL), 'Control Lot' (56777), 'Control ID' (56771), 'Control Lot Expiration' (January 31, 2022 11:59 PM), 'Assay' (ALB, CAL), 'Usage' (Primary), 'Definitions' (CHEMQUAL 56771, CHEMQUAL 56772), and 'Status' (Enable). Buttons for 'Save' and 'Cancel' are at the bottom right of the modal.

Control Name	Control Lot	Control ID	Control Lot Expiration	Assay	Usage	Definitions	Status	Last Order Submitted
CHEMQUAL	56777	56771	January 31, 2022 11:59 PM	ALB, CAL	Primary	CHEMQUAL 56771 CHEMQUAL 56772	Enable	
GluControl	12345	45281	January 31, 2022 11:59 PM	Gluc	Primary	GluControl 12341	Enable	
DRUGQUAL	87659	68541	January 31, 2022 11:59 PM	AFP, CEA, Dig, eE2, hCG, TSH3-UL	Primary	DRUGQUAL 56771 DRUGQUAL 56772	Disabled	
MULTIQUAL	22000	98561	January 31, 2022 11:59 PM	AlbP, AST, CHOL_2, Dgn, K, Na	Primary	MULTIQUAL 78951 MULTIQUAL 78952	Disabled	

	Control Name	Levels	Assay	Control Lot	Control ID	Expiration
> <input type="checkbox"/> ALPQUAL	1,2	Lac, TIBC		32654		
> <input type="checkbox"/> CO2QUAL	1,2,3,4,5	ALB, CAL, CHOL, CREA, Direct Bilirubin		78540		
> <input type="checkbox"/> CHEMQUAL	1,2,3,4,5,6	ALB, CAL		56777		
> <input type="checkbox"/> CRPQUAL	1,2,3,4,5,6	AFP, CEA, Dig, eE2, hCG, TSH3-UL		35968		
> <input type="checkbox"/> DRUGQUAL	1,2,3,4,5,6	AFP, CEA, Dig, eE2, hCG, TSH3-UL		87659		
> <input type="checkbox"/> ENZQUAL	1,2,3,4	Carb, Gent, Tob, Vanc		75130		
> <input type="checkbox"/> GluControl	1,2,3	Gluc		12345		
> <input type="checkbox"/> LYQUICHECK	1,2,3	Na, K, Cl		65214		
> <input type="checkbox"/> MULTIQUAL	1,2,3	AlbP, AST, CHOL_2, Dgn, K, Na		22003		

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Включает:
- Security & Accounts
- Lab Settings
- Calibration и QC Settings
- Patient Results Management
- HIL Configuration
- Alerts, Panels, Ratios
- Test Display Order

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Customer Training

Atellica Solution System Orientation

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Additional Orientation Topics

- **Online Help Orientation**
- **SRS and Software Update Orientation**
- **Emergency Sample Load Orientation**

Customer Training

Review! Customer Training Resources

New CS Product Information Portal:

Atellica Solution

LS > Atellica Solution

Service Knowledge Base Knowledge Base - SKB

Service Parts Catalog

- Atellica IM 1600 Analyzer (11066000)
- Atellica IM 1300 Analyzer (11066001)
- Atellica CI 930 Analyzer (11067000)
- Atellica Direct Load (11068008)
- Atellica Sample Handler Prime (11069001)
- Atellica Sample Handler Additional (11069004)
- Atellica Decapper (11069018)
- Atellica Decapper (11069400)
- Atellica Bone, fat sealer Kit (11069429)

Links

- Document Library (Operator guides, IFUs, customer bulletins, and safety data sheets)
- ConcertMAN (Technical, application, and informatics trainings for service personnel)
- PiConnect (Product overview training for customers and service personnel)
- In-vitro Handover Customer Trainings CB DOC

Download CB-DOC **Open CB-DOC**

Cloud Based Documentation
AP01-000-805.01.01.02 | 2023-09-18 14:18:54 GMT+00:00

in-vitro Handover Customer Trainings

Application

- Laboratory Diagnostics (LD)
 - LD Systems
 - Instructions
 - Application Training Report
- Molecular Diagnostics (MD)
 - Instructions
 - Application Training Report
- Diagnostics IT
 - Instructions
 - Application Training Report
- Atellica Solution
 - Instructions
 - Training Checklist (en)**
 - Application Training Report

Current Product Information Portal:

SIEMENS
Healthineers

Product Information

Atellica Solution (SCI / DL)

Links

- ConcertMAN (Technical, application, and informatics training for service personnel)
- Document Library (Operator guides, IFUs, customer bulletins, and safety data sheets)
- PiConnect (Product overview training for customers and service personnel)

Service Knowledge Base

Service Parts Catalog

CB-DOC

CB-DOC for Download

Customer Training Materials v1.20 and Above

- x 126 to 128 Update Training Materials
- x Advanced Customer Facilitator Guide
- x CI 1980 Bridge Training Materials (for customers experienced on Atellica Solutions)
- x Handover Customer Training Workbook
- x 126 Update Training Materials French
- x Advanced Customer Training Workbook
- x CI 1980 Customer Facilitator Guide
- x Handover Customer Facilitator Guide